


# Tristel Duo™

# ULT

High-level disinfectant foam

FOR NON-INVASIVE  
ULTRASOUND PROBES

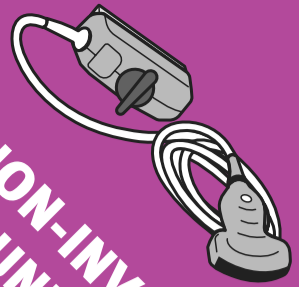

## Notes:

- Do not use if the bottle is damaged.
- Disinfect hands and wear appropriate Personal Protective Equipment (PPE) when handling disinfectants and medical devices.
- If Tristel Duo ULT is used for the first time, depress the pump two to four times to prime it.
- Use Tristel Duo ULT on surfaces that are free from ultrasound gel.
- If the surface is visibly soiled, follow the hospital protocol for cleaning instructions.
- For professional use.
- This wall chart should be used in conjunction with the instructions for use (IFU), product label, Safety Data Sheet (SDS) and device manufacturer's instructions.

## STEP 1



Dispense two aliquots of Tristel Duo onto a dry wipe (Duo Wipes are recommended).

## STEP 2



Use the wipe to spread the foam over the surface of the probe and ensure all areas are covered. Continuous wiping of the surface is not necessary and should be avoided to prevent shedding of the wipe.

**Note:** Discard the used wipe in accordance with local regulations. Do not macerate. Do not re-use.

## STEP 3



Leave surfaces to dry to ensure a minimum contact time of **30 seconds**.

**Note:** Store the medical device in accordance with hospital protocols to prevent damage or recontamination.

**Tristel™**

Created by: Tristel Solutions Limited, Lynx Business Park, Cambs, UK, CB8 7NY  
T +44 (0) 1638 721500 - E mail@tristel.com - W www.tristel.com

Hong Kong & Taiwan: 21st Floor, 168 Electric Road, Fortress Hill, Hong Kong  
T +852 2895 6968 - F +852 2869 4388

For Tristel patent information please visit: <http://www.our-patents.info/tristel>

Copyright © Tristel Solutions  
Mkt-Wal-577-3  
February 2020