
1

**Tristel**™

Pre-Clean Wipe
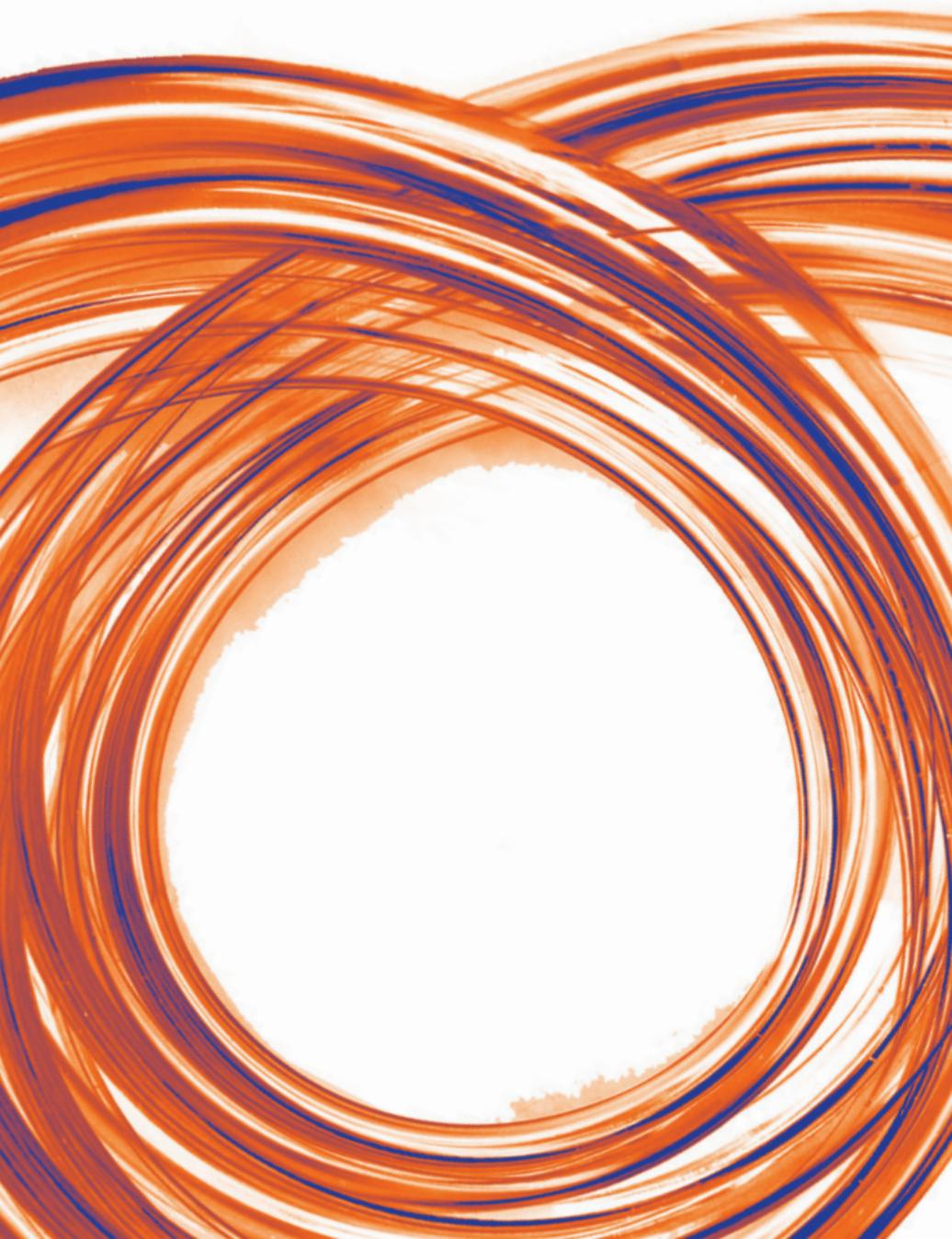
## USER INSTRUCTIONS

- Disinfect hands and wear gloves when handling disinfectants and medical devices.
- Wear appropriate Personal Protective Equipment (PPE).
- Do not use if the Wipe sachet is damaged or the Wipe is discoloured, damaged or dry.
- Keep the empty Wipe sachets to complete the Tristel Quality Audit Trail Record Book at the end of the decontamination procedure.
- Continuous wiping of the surface is not necessary and should be avoided to prevent shedding of the Wipe.
- For professional use only.
- These Instructions for Use (IFU) should be used in conjunction with the product label, Safety Data Sheet (SDS) and device manufacturer's instructions.
- Do not use past the expiry date.

The first step in the decontamination procedure of medical devices is cleaning of the surface to remove soil and organic matter prior to high-level disinfection. The Tristel Pre-Clean Wipe is impregnated with a triple-enzymatic detergent and surfactant.

The Tristel Pre-Clean Wipe is CE Marked as a Class I Medical Device (MDD 93/42 EEC).

### STEP 1

Disinfect hands and wear gloves when handling disinfectants and medical devices.

### STEP 2

Take one Pre-Clean Wipe sachet.

### STEP 3

Remove the Wipe from its sachet and lay it out in the palm of your hand.

### STEP 4

Wipe the surface of the medical device until soil and organic matter have been visibly removed. In case of heavy soiling more than one Wipe may be used.

### STEP 5

Discard the used Wipe and gloves in accordance with local regulations. Do not reuse. Keep the empty Wipe sachet for traceability.

Please complete the Tristel Quality Audit Trail Record Book to keep record of the decontamination procedure.

Manufactured by: Tristel Solutions Limited  
Lynx Business Park, Fordham Road, Snailwell, Cambridgeshire CB8 7NY

T +44 (0)1638 721500 - F +44 (0)1638 721911 - E mail@tristel.com - W www.tristel.com

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Hersteller: Tristel Solutions Ltd
Vertrieb in Deutschland: Tristel GmbH, Karl-Marx-Allee 90A, 10243 Berlin, Deutschland
Vertrieb in der Schweiz: Tristel AG, Sandgrube 29, CH-9050, Appenzell
Tristel Solutions Limited, Lynx Business Park, Fordham Road, Snailwell, Cambridgeshire, CB8 7NY
**Australian sponsor:** Tristel Pty Ltd. 40/328 Reserve Road, Cheltenham, VIC 3192, Australia