
3

Tristel™

Rinse Wipe

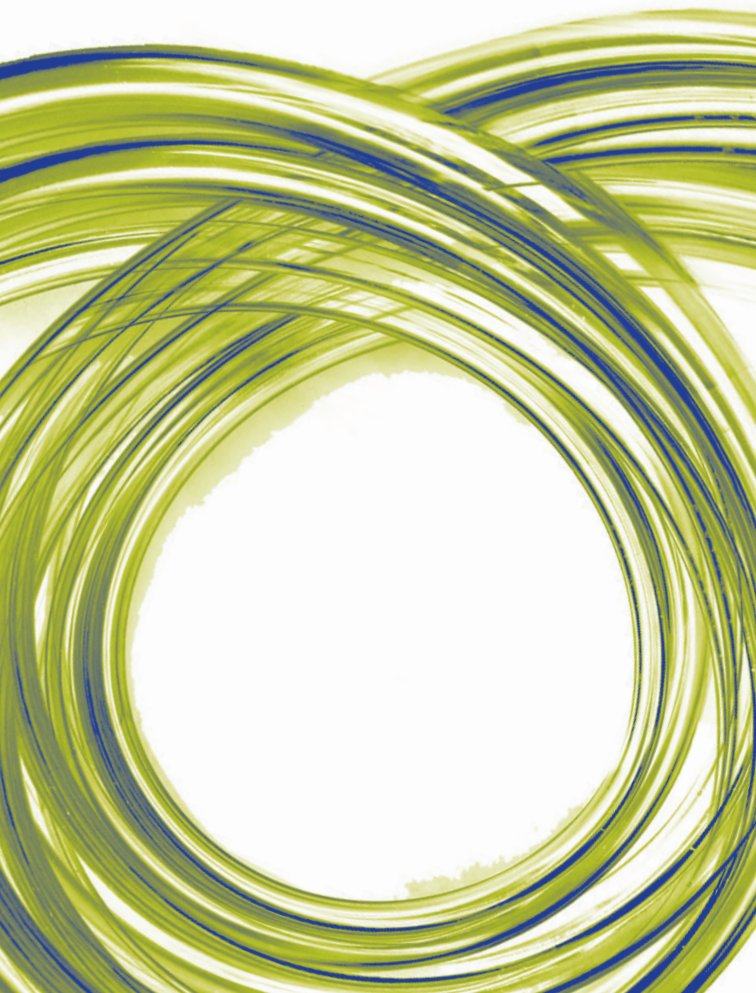USER INSTRUCTIONS

- Disinfect hands and wear gloves when handling disinfectants and medical devices.
- Wear appropriate Personal Protective Equipment (PPE).
- Do not use if the Wipe sachet is damaged or the Wipe is discoloured, damaged or dry.
- Keep the empty Wipe sachets to complete the Tristel Quality Audit Trail Record Book at the end of the decontamination procedure.
- Continuous wiping of the surface is not necessary and should be avoided to prevent shedding of the Wipe.
- For professional use only.
- These Instructions for Use (IFU) should be used in conjunction with the product label, Safety Data Sheet (SDS) and device manufacturer's instructions.
- Do not use past the expiry date.

The third and final step in the decontamination procedure is rinsing of the medical device. The Rinse Wipe is impregnated with de-ionised water and a low-level of antioxidant which removes chemical residues from a surface. Each Rinse Wipe sachet is packed and then sterilised by gamma-irradiation.

The Rinse Wipe is CE Marked as a Class I Sterile Device (MDD 93/42 EEC).

**STEP 1**

Take one Rinse Wipe sachet.

STEP 2

Remove the Wipe from its sachet and lay it out in the palm of your hand.

STEP 3

Wipe the surface of the device that has been decontaminated to remove excess foam.

STEP 4

Discard the used Wipe and gloves in accordance with local regulations. Do not reuse. Keep the empty Wipe sachet for traceability.

Note: Upon completion of the decontamination cycle the device should be left to air dry. Store the device in accordance with hospital protocols to prevent damage or recontamination.

Please complete the Tristel Quality Audit Trail Record Book to keep record of the decontamination procedure.

Manufactured by: Tristel Solutions Limited
Lynx Business Park, Fordham Road, Snailwell, Cambridgeshire CB8 7NY
T +44 (0)1638 721500 - F +44 (0)1638 721911 - E mail@tristel.com - W www.tristel.com

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The Tristel Rinse Wipes is classified under the Australian Register of Therapeutic Goods as a Class Is medical device, AUST R number: 173046.

For Tristel patent information please visit www.our-patents.info/tristel

Contact Tristel, your local distributor or visit www.tristel.com for supporting documents such as safety data sheets, microbiological test data and reports.

Tristel NV, Smallandlaan 14B,
2660 Antwerp, Belgium

RIN/098/4
December 2020

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Hersteller: Tristel Solutions Ltd
Vertrieb in Deutschland: Tristel GmbH, Karl-Marx-Allee 90A, 10243 Berlin, Deutschland
Vertrieb in der Schweiz: Tristel AG, Sandgrube 29, CH-9050, Appenzell
Tristel Solutions Limited; Lynx Business Park, Fordham Road, Snailwell, Cambridgeshire, CB8 7NY
Australian sponsor: Tristel Pty Ltd. 40/328 Reserve Road, Cheltenham, VIC 3192, Australia