


THE ONLY HIGH-LEVEL DISINFECTANT PROVEN EFFECTIVE AGAINST HPV IN 30 SECONDS, DESIGNED TO TREAT THE ENTIRE ULTRASOUND STATION

There are over 200 different Human Papillomavirus (HPV) types. Types 16 and 18 are known to be responsible for approximately 70% of cervical cancer cases. A new scientific study by Meyers, et al. (2020) demonstrates that Tristel Duo ULT possesses efficacy against infectious HPV types 16 and 18 on a transvaginal ultrasound probe in 30 seconds.

Tristel Duo ULT provides high-level disinfection for endocavity ultrasound probes, skin surface ultrasound probes (including those used during invasive procedures), probe holders, cords, plugs, keyboards, monitors and even mobile medical devices that combine an ultrasound probe with a smart device.

TRISTEL DUO ULT IS EASY TO USE. SIMPLY APPLY, WIPE AND WAIT!

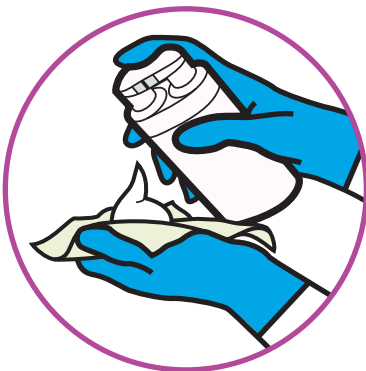

DISPENSE TWO ALIQUOTS OF
TRISTEL DUO ULT ONTO A DRY WIPE
(DUO WIPE IS RECOMMENDED).



USE THE WIPE TO SPREAD THE FOAM
OVER THE SURFACE OF THE ULTRASOUND
PROBE AND ASSOCIATED PARTS.
ENSURE ALL AREAS ARE COVERED.



LEAVE THE SURFACE TO DRY TO
ENSURE A MINIMUM CONTACT TIME
OF 30 SECONDS.

Tristel DuoTM
ULT



EFFECTIVE

Tristel Duo ULT is sporidical, mycobactericidal, virucidal, fungicidal and bactericidal with a contact time of 30 seconds. The product has been extensively tested according to European Standard suspension tests and a standardised methodology that involves the inoculation of surfaces with the test organism. Tristel's chlorine dioxide is effective against microorganisms of concern such as:

- *Bacillus subtilis*
- *Mycobacterium terrae*
- **Human Papillomavirus (HPV)**
- Hepatitis B Virus (HBV)
- Human Immunodeficiency Virus (HIV)
- Herpes Simplex Virus Type 1
- *Aspergillus brasiliensis*
- *Candida albicans*
- *Escherichia coli*
- *Neisseria gonorrhoeae*
- *Staphylococcus aureus*
- *Gardnerella vaginalis*
- *Streptococcus agalactiae*
- Carbapenem-resistant Enterobacteriaceae (CRE)
- *Klebsiella pneumoniae*
- Vancomycin-resistant Enterococci (VRE)
- *Enterococcus faecium*
- Methicillin-resistant *Staphylococcus aureus* (MRSA)



SIMPLE

Training and certification are an essential part of Tristel's customer service. Our Online Training Portal provides Tristel Duo ULT user training at a time and place convenient to you.



ECONOMICAL

Tristel Duo ULT does not require investment in expensive equipment or service contracts. Two aliquots of foam and one Duo Wipe is all it takes to achieve disinfection to the highest level.



MOBILE

Tristel Duo ULT can be used wherever required, without having to worry about the limitations that automated systems bring. Tristel Duo ULT is also suitable for use in remote locations where electricity and water are not readily available.




COMPATIBLE

Tristel Duo ULT has been tested and confirmed to be compatible by major manufacturers, including:

- Alpinion Medical
- BD (Bard Access)
- BK Medical
- Butterfly Network
- Canon Medical Systems
- Carestream
- Fujifilm (Hitachi)
- FUJIFILM SonoSite
- GE Healthcare
- Healcerion
- MCube
- Mindray
- MobileODT
- Philips
- Quantel Medical
- Samsung Healthcare
- Siemens Healthineers
- Sonoscape
- Supersonic Imagine
- Verathon

APPLICATIONS

Tristel Duo ULT is the only high-level disinfectant that treats all parts of an ultrasound probe the same. It is the only disinfectant that can travel alongside mobile medical devices as an easy-to-use high-level disinfection method.

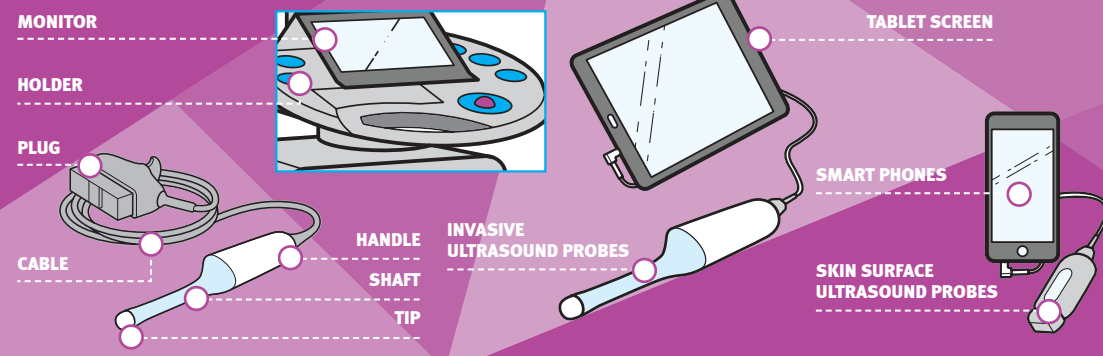

Duo Wipes are designed specifically for the application of Tristel Duo ULT.

Duo Wipes are made of 100% non-woven, low-linting polypropylene (18gsm) and pre-cut to prevent shedding when dispensed from the tub.

The use of Duo Wipes allows the effective release of the disinfectant onto the surface.



DUO WIPES

- 200 wipes per tub / 6 tubs per box

TRISTEL DUO ULT

- 125ml Tristel Base Solution (citric acid) + 125ml Tristel Activator Solution (sodium chlorite)

When mixed upon pressing the foam pump, Tristel's proprietary chlorine dioxide chemistry is generated.

- 310 aliquots per bottle / 6 bottles per box

DUO DOK

- Designed to keep Tristel Duo ULT to hand
- 5 per pack

Tristel

Medical Device
Decontamination

UK Manufacturer: Tristel Solutions Limited, Lynx Business Park, Cambs, UK, CB8 7NY
T +44 (0) 1638 721500 - E mail@tristel.com - W www.tristel.com

Malaysia: Tristel Malaysia Sdn. Bhd, A-25-3A Pinnacle PJ, Lorong Utara C, 46200 Petaling Jaya, Selangor Darul Ehsan, Malaysia
T +603 79317791 - F +603 79319971 - E mail-my@tristel.com - W www.tristel.com/my-en

For Tristel patent information please visit: <http://www.our-patents.info/tristel>

Copyright © Tristel Solutions
Mkt Bro 1504-4
10/NOV/2022

Tristel
WE HAVE CHEMISTRY.