


INTERMEDIATE-LEVEL DISINFECTION FOR NON-INVASIVE ULTRASOUND PROBES

NOTES

- Do not use if the bottle is damaged.
- Disinfect hands and wear suitable gloves and Personal Protective Equipment (PPE) when handling disinfectants and medical devices.
- When Tristel DUO NCU is being used for the first time, depress the pump two to four times to prime it.
- If the surface is visibly soiled, follow the hospital protocol for cleaning instructions.
- For healthcare professional use only.
- This wall chart should be used in conjunction with the instructions for use (IFU), product label, Safety Data Sheet (SDS) and device manufacturer's instructions.


DIRECTIONS FOR USE

1

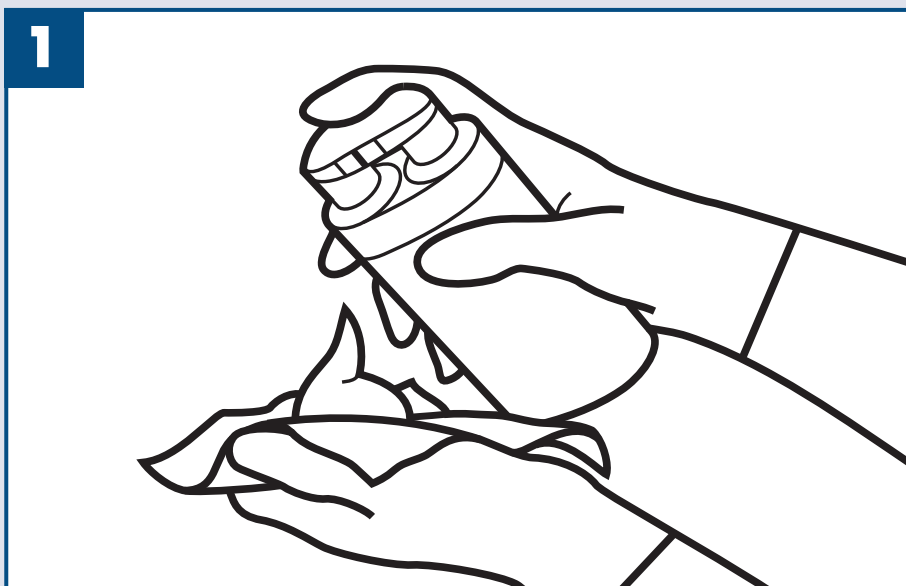


Depress the pump to dispense Tristel DUO NCU onto a dry wipe (DUO WIPES are recommended) or directly onto the surface.

2

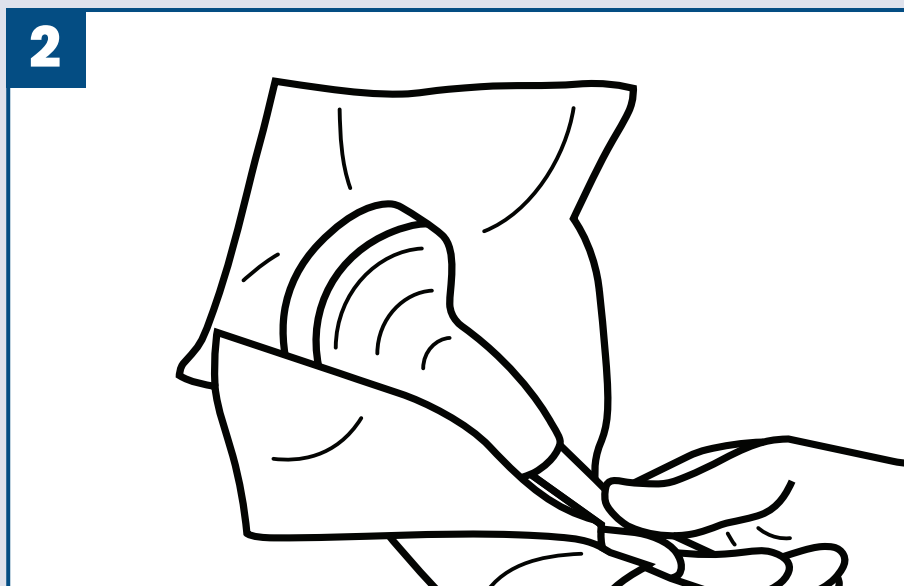


Use the wipe to spread the foam over the surface and ensure all areas are wet.

3

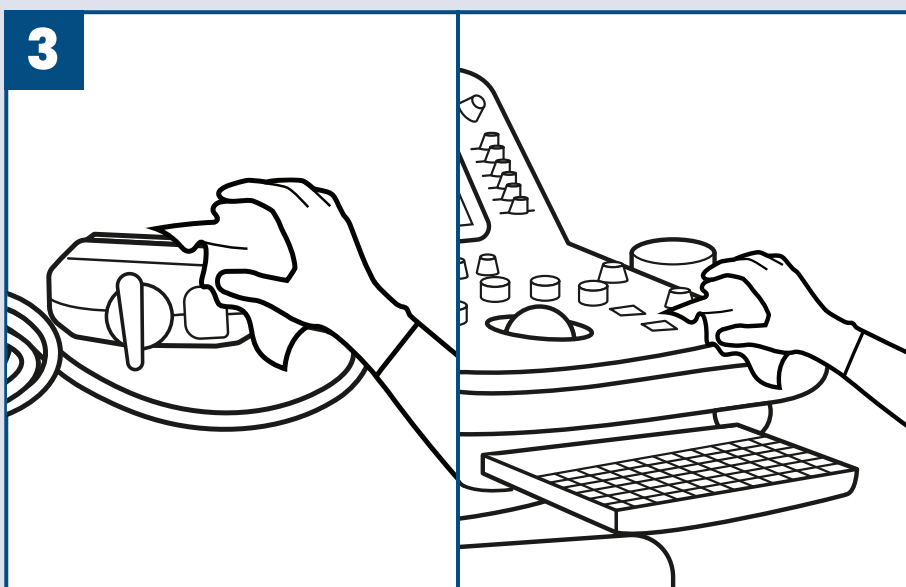

Use Tristel DUO NCU to disinfect surfaces that can come into contact with the probe, such as holders, cables and plugs.

4

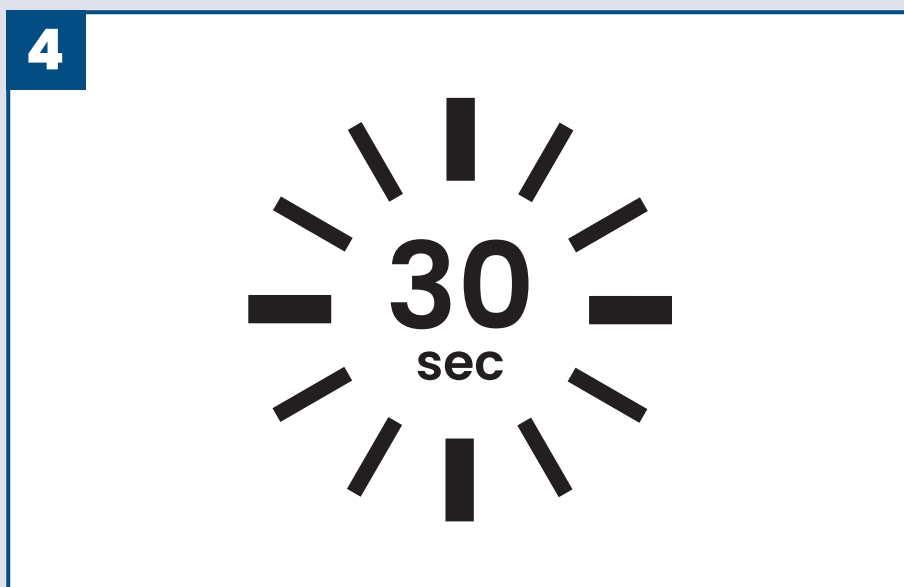

Allow the surface to remain wet for 30 seconds. Leave to air dry.

Note: Discard the used wipe to clinical waste. Do not macerate. Do not re-use.