


Tristel Duo™

OPH

FOR REUSABLE
TONOMETER PRISMS

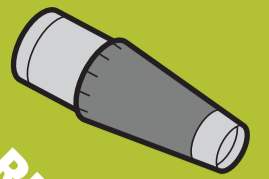

Notes:


- Do not use if the bottle is damaged.
- Disinfect hands and wear appropriate Personal Protective Equipment (PPE) when handling disinfectants and medical devices.
- Remove and discard transport locks that stop the pump from being depressed in transit.
- If Tristel Duo OPH is used for the first time, depress the pump two to four times to prime it.
- If the surface is visibly soiled, follow the hospital protocol for cleaning instructions.
- For professional use.
- This wall chart should be used in conjunction with the instructions for use (IFU), product label, Safety Data Sheet (SDS) and device manufacturer's instructions.

STEP 1



Dispense two aliquots of Tristel Duo onto a dry wipe (Duo Wipes are recommended) or directly onto the prism.

STEP 2

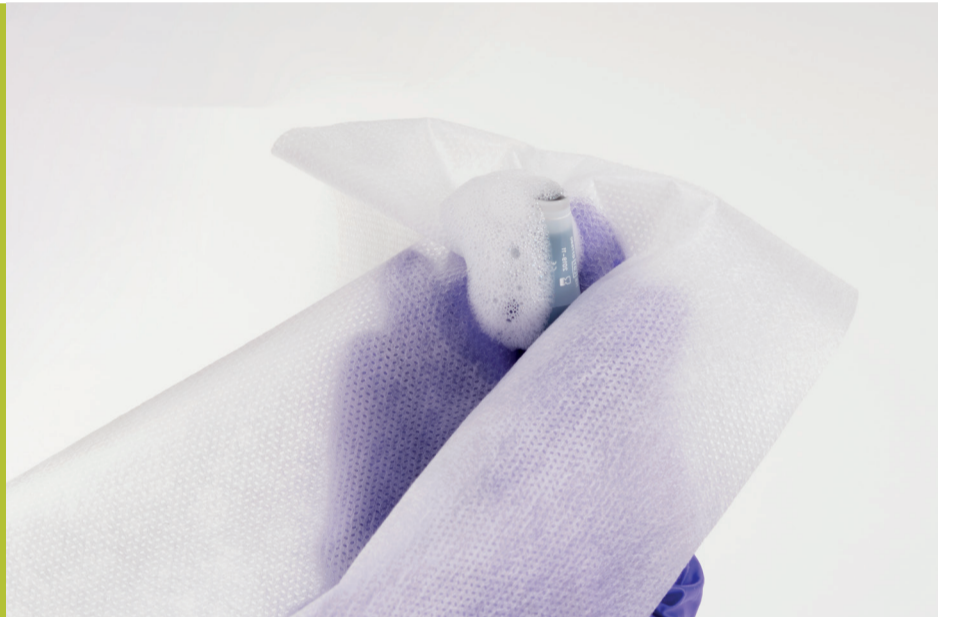

Use the wipe to spread the foam over the surface of the lens and ensure all areas are covered. Continuous wiping of the surface is not necessary and should be avoided to prevent shedding of the wipe.

Note: Discard the used wipe in accordance with local regulations. Do not macerate. Do not re-use.

STEP 3



Leave the surface to dry to ensure a minimum contact time of **two minutes**.

STEP 4



Rinse the lens with water of appropriate quality.

Note: Store the medical device in accordance with hospital protocols to prevent damage or recontamination.

Tristel™

Created by: Tristel Solutions Limited, Lynx Business Park, Cambs, UK, CB8 7NY
T +44 (0) 1638 721500 - E mail@tristel.com - W www.tristel.com

For Tristel patent information please visit: <http://www.our-patents.info/tristel>

Australian sponsor: Tristel Pty Ltd, 40/328 Reserve Road, Cheltenham, VIC 3192
T 1300 680 898 - F +61 (0)3 9533 6193 - E mail-au@tristel.com

New Zealand sponsor: Tristel New Zealand Limited, P.O Box 15451, Tauranga, New Zealand
T +64 (0) 7 5771560 - E info@tristel.co.nz

Copyright © Tristel Solutions
Mkt-Wal-350-2
March 2019