

# Tristel Trio™

## Wipes System


A complete solution for high-level disinfection of non-lumened invasive medical devices.

 Refer to the User Guide for full instructions.

### Notes:

- This wall chart should be used in conjunction with the instructions for use (IFU), product label, Safety Data Sheet (SDS) and device manufacturer's instructions.
- Always follow hospital protocols and device manufacturers' instructions.
- Disinfect hands and wear appropriate PPE when handling disinfectants and medical devices.
- An initial pre-clean should be performed immediately after medical device use as per hospital protocol.
- Take care not to rub or bend the instrument and pay particular attention to edges, ridges, indentations and areas where different materials connect.
- Prepare for the decontamination process by taking one of each Wipe and the Activator Foam from the box.
- Do not use past the expiry date. For the expiry date and LOT number please see the surface of the carton.
- Do not use if the Wipe sachet is damaged, or the Wipe is discoloured, damaged or dry.
- Use more Pre-Clean Wipes if required to ensure the instrument is visibly clean.
- Keep the Wipe sachets to complete the Tristel Traceability book at the end of the decontamination procedure.
- Discard the used Wipe and gloves in accordance with local regulations.
- For professional use only.
- Establish a dirty and clean area. Adhere to local and hospital decontamination protocol when establishing a clean and dirty area.

FOR FLEXIBLE  
NON-LUMENED CYSTOSCOPES

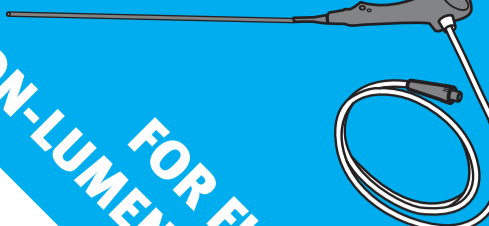

**STEP 1 : CLEANING**



Remove the Pre-Clean Wipe from its sachet and unfold it in the palm of your hand.

Wipe the probe to remove all visible soil.


Place the instrument in the dirty area.

Discard the used Wipe and gloves in accordance with local regulations.

Disinfect your hands and put on new gloves.



**STEP 2a + 2b : ACTIVATING**




Remove the Sporidical Wipe from its sachet and unfold it in the palm of your hand.

Take the Activator Foam bottle and apply two doses of foam onto the Sporidical Wipe.


Fold the Sporidical Wipe in on itself and scrunch together **15 seconds** to generate chlorine dioxide.

Ensure that the Wipe is completely covered in Foam.

Use the activated Wipe immediately.



**STEP 2 : HIGH-LEVEL DISINFECTING**



Wipe the surface of the probe in one movement to cover it with foam, ensuring all areas come into contact with the Wipe.

Discard the used Wipe in accordance with local regulations.

Observe a minimum contact time of **30 seconds**.

Place the instrument in the clean area.



**STEP 3 : RINSING**



Remove the Rinse Wipe from its sachet and unfold it in the palm of your hand.

Wipe the probe to remove excess foam.

Discard the used Wipe, apron and gloves in accordance with local regulations. Disinfect your hands.



**STEP 4 : TRACEABILITY**



The Tristel Trio Wipes System can include either paper-based or digital traceability.

- The Tristel Quality Audit Trail Record Book can be used to manually document the Tristel Trio Wipes System decontamination procedure.
- Tristel 3T is a digital traceability system comprising an online Portal for data management and reporting, and an App to record Tristel Trio Wipes System decontamination procedures. The 3T App features optional training videos for increased compliance.

Both traceability systems include guidance on how to complete the record process correctly.