

# Tristel Duo™

# OPH

FOR DIAGNOSTIC  
CONTACT LENSES



High-level disinfectant foam

## Notes:

- Do not use if the bottle is damaged.
- Disinfect hands and wear appropriate Personal Protective Equipment (PPE) when handling disinfectants and medical devices.
- If Tristel Duo OPH is used for the first time, depress the pump two to four times to prime it.
- If the surface is visibly soiled, follow the hospital protocol for cleaning instructions.
- For professional use.
- This wall chart should be used in conjunction with the instructions for use (IFU), product label, Safety Data Sheet (SDS) and device manufacturer's instructions.

### STEP 1



Dispense two aliquots of Tristel Duo onto a dry wipe (Duo Wipes are recommended) or directly onto the lens.

### STEP 2



Use the wipe to spread the foam over the surface of the lens and ensure all areas are covered. Continuous wiping of the surface is not necessary and should be avoided to prevent shedding of the wipe.


**Note:** Discard the used wipe in accordance with local regulations. Do not macerate. Do not re-use.

### STEP 3



Leave the surface to dry to ensure a minimum contact time of **30 seconds**.

### STEP 4

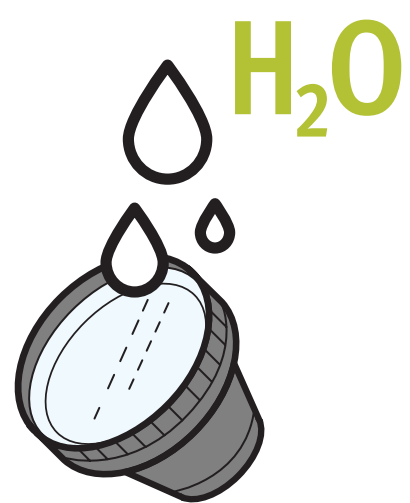

Rinse the lens with water of appropriate quality.

**Note:** Store the medical device in accordance with hospital protocols to prevent damage or recontamination.

**Tristel™**

**Created by:** Tristel Solutions Limited, Lynx Business Park, Cambs, UK, CB8 7NY  
T +44 (0) 1638 721500 - E [mail@tristel.com](mailto:mail@tristel.com) - W [www.tristel.com](http://www.tristel.com)

For Tristel patent information please visit: <http://www.our-patents.info/tristel>  
Copyright © Tristel Solutions - Mkt-Wal-004-9 - 14/APR/2022

**Hong Kong & Taiwan:** Tristel Asia Limited, 21st Floor, 168 Electric Road, Fortress Hill, Hong Kong  
T +852 2895 6968 - F +852 2869 4388 - E [customerservicehk@tristel.com](mailto:customerservicehk@tristel.com)

**Singapore:** Tristel Pte Ltd, 2 Venture Drive, Vision Exchange #14-04, Singapore 608526  
T +65 6518 9918 - E [mail-sg@tristel.com](mailto:mail-sg@tristel.com)

**Malaysia:** Tristel Malaysia Sdn. Bhd, A-25-3A Pinnacle PJ, Lorong Utara C, 46200 Petaling Jaya, Selangor Darul Ehsan, Malaysia  
T +603 79317791 - F +603 79319971 - E [mail-my@tristel.com](mailto:mail-my@tristel.com)